

# InnoCyto Inc.

15375 Barranca Pkwy, Suite I-103 Irvine, CA 92618

## **Technical Data Sheet**

iF647 Anti-human CD34

Catalog Number: 107004, 107005

Size: 25 tests, 100 tests

Target Name: CD34, Gp105-120, My10

Regulatory Status: RUO

### **Product Details**

Clone: 034AB
Application: FC
Reactivity: Human
Format: iF647
Isotype: Rabbit IgG

Antibody Type: Monoclonal

Formulation: Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and 0.2% (w/v)

BSA

**Protein Concentration:** Supplied at a lot-specific concentration.

Storage&Handling: The antibody solution should be stored undiluted between 2°C and 8°C, and

protected from prolonged exposure to light. Do not freeze.

Recommended Usage: For flow cytometric staining, it is recommended to use 5 uL of this reagent per 0.5-1.0 million cells in a 100  $\mu$ L volume. Optimal reagent performance should be determined by

titration for each specific application. **Excitation Laser:** Red Laser (633 nm)

Release Date: Dec-25

## **Background Information**


CD34 (gp105–120) is a type I monomeric sialomucin-like glycophosphoprotein with a molecular weight of 105–120 kDa, expressed on hematopoietic stem and progenitor cells, bone marrow stromal cells, capillary endothelial cells, embryonic fibroblasts, and some neural tissues. CD34 serves as a widely used stem cell marker and mediates cell adhesion and lymphocyte homing through interactions with L-selectin and E-selectin ligands. CD34 expression is dynamic and reversible, reflecting specific states of hematopoietic development, including proliferation, differentiation, and altered adhesion properties. Multiple epitopes have been identified based on differential enzymatic sensitivity. Clinically, CD34 is expressed in ~40% of acute myeloid leukemias and 65% of pre-B acute lymphoblastic leukemias, but only 1–5% of acute T-lymphoid leukemias. It is also used to assess stem cell counts for transplantation and can provide prognostic information in certain cancers.

### **Product Data**



# InnoCyto Inc.

15375 Barranca Pkwy, Suite I-103 Irvine, CA 92618

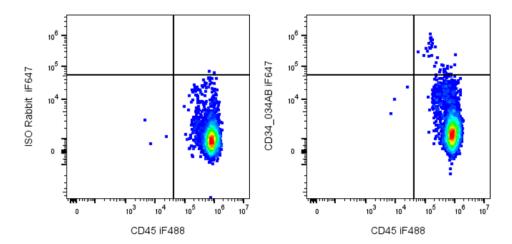

Human peripheral blood leukocytes stained with iF488 Anti-human CD45 and either FITC Anti-Human CD34 clone 034AB (right panel) or an isotype control (left panel).